
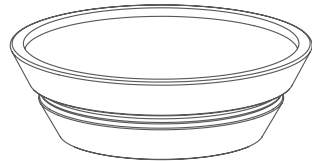


Manual Petri Dish Spreader

DS-M

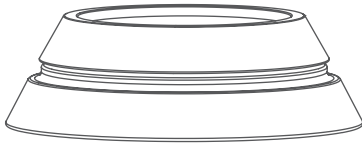

1.Introduction

This product is designed for manual spreading of agar plate in petri dish, which is mainly used in the fields of molecular biology and microbiology.

2.Features

- Manual operation, easy to control;
- Antiskid design, easy to handle;
- Double sides usable, two work spaces in one device;
- Compatible with different sizes of petri dish.

3.Technical Data

Cat.No	DS-M
Petri Dish Size (Φ)	90-150mm
Rotate Mode	Manual
Work Plate	2 Sides
Dimensions (Φ*H)	160 * 50 mm
Weight	~2.4kg

4.Operation

Place culture plate on the suitable work space, add sample and rotate the work platform to spread evenly with the L-spreader.

5.Warranty

One year limited warranty. It is not suitable for damage caused by improper use.



Beijing Labgic Technology Co., Ltd.

Add:No.9 Yumin Street, Area B of the Airport Industrial Zone, Shunyi District, Beijing 101318 China



Toll Free:400-600-4213

Website:www.labgic.com



手动平板涂布器

DS-M

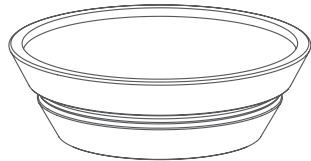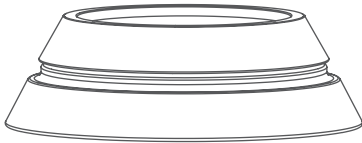

1. 产品简介

本产品用于培养皿内琼脂平板涂布操作，主要应用分子生物学和微生物学等领域。

2. 产品特点

- 手动操作，易控制；
- 简易防滑设计，易操作；
- 双面可用，适用不同规格培养皿。

3. 技术参数

货号	DS-M
培养皿直径 (Φ)	90-150mm
操作模式	手动
工作台	2个
尺寸 (Φ* <i>h</i>)	160 * 50 mm
重量	~2.4kg

4. 使用说明

将培养平板放置到适合的工作台上，滴加样品，旋转工作平台，用涂布棒完成涂布。

5. 保修说明

质保一年，因人为使用不当造成的损坏，不属于质保范围。



北京兰杰柯科技有限公司

地址：北京市顺义区空港工业园区B区裕民大街9号1幢206室

电话：400-600-4213

网址：www.labgic.com

